

The Eliminator Electro-pulsing Blood Cleanser



The Eliminator Blood Pulsing device is precisely based on the original design presented as a free, self-help healing technology to the world by Dr Robert C. Beck in the mid 1990's as the first part of the "Beck Protocol" (consisting of four elements in total). Depending on your state of health you may only need to incorporate this device as a preventative solution and can expect protection against virus and bacteria infections. As a health maintenance tool you may opt for even a more intensified protocol when combined with Colloidal Silver intake using the Colloidal Silver generator with current control and auto-shut-down functionality. The complete Beck Protocol consist of four elements which in totality is targeted as a further curing protocol for people with serious ailments and has shown excellent results for a variety of health-related conditions including HIV/Aids, and other viral or bacterial infections. These treatments when used within recommendation limits have ben found toi have little or insignificant side-effects. It migh be of interest to also investigate and combine these devices/treatments with another very safe and effective substance namely DMSO as an even more intense treatment protocol for a wide variety of serious health conditions.

The original design of the Blood-Cleansing Pulsar revealed by Dr Robert Beck is published here: <https://www.bobbeck.com>

I encourage all users of this device to visit the Beck Protocol webpage, watch the videos, read the testimonials and download the PDF instruction booklet on the complete original Beck Protocol here: <https://www.bobbeck.com/pdfs/guide-to-using-the-beck-protocol.pdf>

Watch this interesting video on YouTube (while available): <https://youtu.be/xMlXd9hLpcY> (or search for "Dr Robert Bob Beck Lecture - Suppressed Medical Discovery To Treat All Disease" the Granada Forum).

You can construct your own unit following the information from the link above. However, in addition to the basic circuit and specifications originally released by Dr Beck, the Eliminator Electro-Pulsing Blood Cleanser offers these added functions:

1. The Eliminator can be used with different power sources:
A 240V AC / 5V DC adaptor that is supplied with the unit as an alternative power source. The unit can also be powered by either a 9V (PP3) battery or two 1.5V (AA size) batteries. Both of these options will fit in the battery compartment.
2. Detection circuitry has been added to visually confirm connection to the body and to measure the peak-pulse current flow. The unit provides visual indication of the 3 -7ma range as specified in the Beck protocol.
3. The eliminator can work from any power supply (DC adaptor or batteries) from 3V to 24V with automatic internal voltage step-up circuitry regulating the voltage at a constant 30V and there is no need for the three series 9V batteries as per original Beck design.
4. The colloidal silver generator (optionally available) offers 1mA current control which is essential to ensure the smallest particle size silver particles for best body cell absorption. It also offers an auto-shutdown feature calibrated at ~ 5 -10ppm to ensure well below toxic levels are maintained.

Items included with the Eliminator Electro-pulser:

1. The Eliminator Electro-Pulsing Handheld Unit
2. Two elastic wrist straps with fixed electrode sleeves
3. One 240VAC / 5V DC mains adaptor



4. Optional (Not included) : Two 1.5V batteries (size AA) or 9V battery

As a non-medical professional I am not allowed to make any definite claims, offer guarantees and advise on medical conditions in regards to your specific health outcome, but

with this well-tested technology, myself, family and friends can testify to the effectiveness of this device. You may improve your health and sometimes experience the immediate effects when treating bodily pains and other afflictions. Many users reported that they do not require vaccinations, medicines for colds and a myriad of other viral or bacterial infections anymore. Based on the research by Dr Beck and other medical professionals who investigated electro-pulsing technology, this may just as well offer the broadest spectrum of protection and cure against the more than 640 viruses known to human kind and all harmful bacteria that invaded our bodies from childhood, using it as a host to feed and multiply. Remember that prevention is always the best cure!

Everyone using the Eliminator Electro-blood-pulsing is expected to experience some health benefits very quickly, even if not knowingly afflicted.

PLEASE NOTE:

Most users will experience side-effects related to detoxification of the body during the first few days after starting to use this device. Reduce pulse time to 10 minutes daily or pulse every second day if necessary until the side effects have cleared up, drink lots of water. Side effects due to detoxification is normal and a good sign that your immune system is starting to expel viruses, bacteria and other pathogens from your body. If kidneys or liver becomes painful, reduce the pulsing time. Other slight symptoms may include a slight headache and/or running nose. These symptom will normally clear up after a few days. You may then increase the pulse time until an hour per day is achieved. After starting out with this treatment, use on a daily basis for at least 30 days to ensure spanning the incubation time of bacteria and viruses in the blood.

If you do not experience better health and vitality after three months of daily use, just return the item in good form for a refund (less delivery cost and any physical damage or shortages).

Why Micro pulse? For Your Health & Quality of Life.

Reduces pathogens and parasites - gentle micro currents of electricity essentially cleanse the blood of bacteria, viruses, fungi, and other pathogens...leading the way to health and vitality.

Increase energy - adenosine triphosphate (ATP) is a key factor in energy at the cellular level. Electricity boosts ATP levels which, in turn increases cellular energy.

Strengthen the body's immune system - micro currents support the body's electrical systems, strengthening the immune system and accelerating the healing process.

Help support cellular function - both red and white blood cells respond to micro currents by operating more efficiently and effectively. Restored voltage helps cells to function at optimal, healthy, normal levels.

Far-reaching benefits - micro current therapy impacts the entire body at a cellular level, using circulation as a delivery system for supplying electrical nutrition to tissues and cells.

Calming - micro currents calm the sympathetic nervous system helping to reduce internal stress levels. It is also reported to clear the mind and improve concentration.

How to use the Eliminator

The Eliminator unit is housed in an ergonomic handheld translucent blue plastic casing that can be carried around, strapped to the arm or inserted into a pocket or pouch. The 1.2m wires with stainless steel electrodes can then be strapped to the forearm/s, ankles or combination of the two by means of the stretchable wrist bands. There are two methods by which the unit can be utilized:

1. By placing the electrodes on main arteries on the wrist/s or ankle/s to cleanse the blood from viruses, bacteria and relief from many other ailments.
2. Localised treatment of any part of the body by positioning the electrodes to both sides of the affected area or on both sides of the body for through-body pulsing of deeper internal affected areas. Current flow can be directed directly through the affected area or organ where possible. (Effective to treat all sorts of bodily infections, pains, arthritis, acne, blisters or sores caused by herpes, dermatitis or stomach pains etc.) Many of these conditions will however also disappear by the normal blood pulsing method on the wrist/s. For localised treatment, make your own smaller electrode sleeves, or just wrap a strip of soft material in two to three layers around each electrode and wet thoroughly.

Although the current level is similar to that of a heart pacer, for theoretical safety reasons we do not recommend the use of this unit for pregnant women, patients with a heart pacer or pulsing directly through the chest and through the heart or directly through the brain, although no-one has reported negative side any effects in this regard.

What to do:

Do blood electrification for a minimum of 1-2 hours per day for a minimum of four weeks to ensure the incubation periods of bacteria and/or viruses have been fully spanned without skipping a day. If a day was skipped, use on a daily basis for another four weeks. Use on a daily basis or occasionally as you prefer after the first four weeks as a preventative measure. I personally recommend using the device on a daily basis for 3 months and you will notice that apart from the absence of colds and flu, you may also experience heightened vitality and energy. Use on a daily basis during cold and flu seasons or epidemics.

Apply normal tap water or salt-water to thoroughly wet the small electrode sleeves on the wrist straps. Position the electrodes over the Ulnar and Radial arteries on opposite insides of same wrist. (The developer of the Eliminator himself is using the electrodes attached to both wrists and have done this every day for many months, and still do, for more effective current flow through the arteries. Pulsing from both wrists does not allow for part of the current to flow (and be lost) through the skin tissue between electrodes. If you decide to use this method, one hour pulsing per day is sufficient due to more effective current flow through the blood.)

In order not to be frowned upon by the health authorities or critics looking for a reason to criticize this protocol, consider the one arm positioning of the electrodes as official, safest and recommended way as presented by Dr Beck.

Velcro® and elastic straps hold the electrodes in place. This should not interfere with other activities. As your blood circulates normally, enough will be flowing along this path in the forearm (or through both arms) until most blood in your body is eventually treated. About 3 to 6mA (milli-ampere) is required at skin contact level to overcome resistive losses through body tissue before current reaches blood within the 50-100uA range.

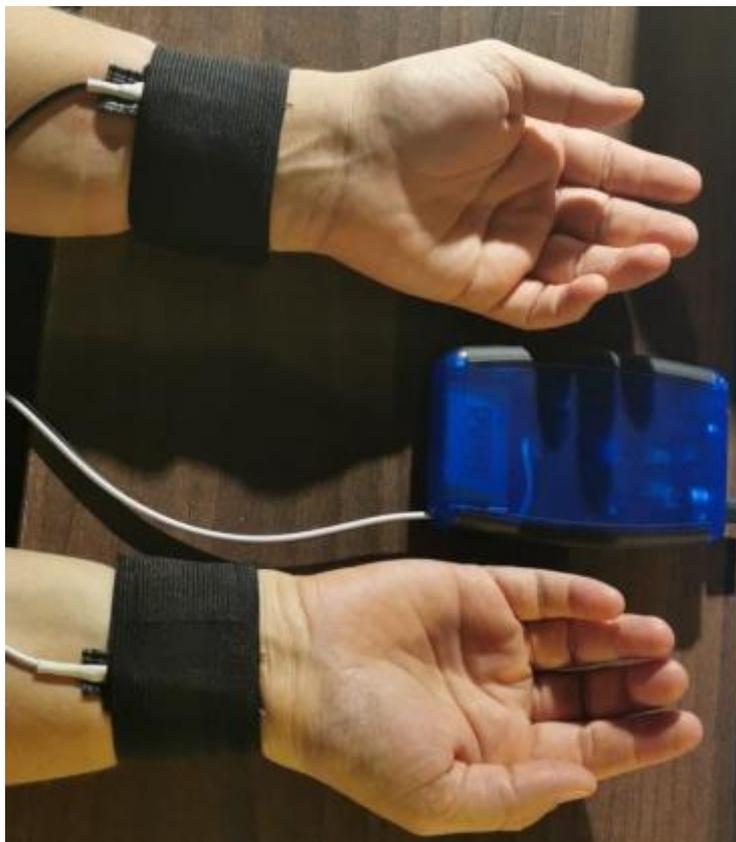



Fig.1 Using the two arms method

To use on one wrist, use only one wrist strap. Thoroughly wet the sleeve with normal tap water and insert one electrode into the fixed electrode sleeve and position on the artery near the thumb where you can detect pulse. Use any soft material strip, wrap around the other electrode and wet thoroughly, push the loose sleeve and electrode underneath the strap on the other main artery on the same wrist at the pinkie's side where you can feel a pulse, see images below for the location of the main arteries on the arm where you can detect a pulse.

If a slight electrical impulse could not be felt during the session, add a little bit of pure salt (not table salt) to the water used to wet the electrodes. If you cannot feel a pulse, watch the LED and adjust to within the range where the internal LED will flicker and up to maximum setting.

I have learned from professional medical documents that the main arteries from the arms are connected above the heart and current flow will not pass through the heart. Electrical current which is only in the micro-amperes range in this case, will always follow the shortest conductive path between electrodes. I have not experienced any side effects from using the electrodes on both arms on a daily basis over many months of using the device.

The Front Panel

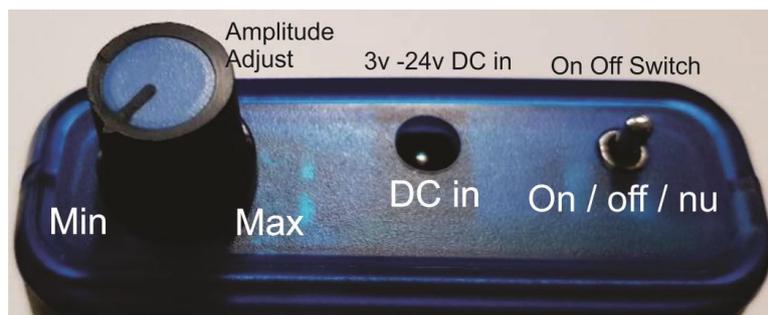

Amplitude adjust – Turn the potentiometer anti-clockwise to Minimum setting before starting your pulsing session, gradually turn clockwise until the LED inside the box starts to flicker, adjust clockwise until a detectable impulse can be felt, or to maximum while no pulse can be felt. The LED should flicker to indicate proper contact and current flow.

Switch

On – the LED will be lit. Move the toggle switch to the left “on” position to start your pulse session).

Off – Middle position is the “off” position, switch the unit off after use to save battery power. Note: the LED will always be on if a DC adaptor is connected.

nu – (Normally not used), This position is normally not used but may be used while connected to the mains adaptor to re-charge a rechargeable 9V battery of good quality. Charge for 12-15 hours. (Use only with a rechargeable 9V battery installed) See section below for charging batteries. Best results with battery use will be achieved using two full capacity Duracell 1.5V batteries or 1.5V rechargeable batteries which should be separately charged with a dedicated charger.

Follow these steps to apply the Eliminator Electro-pulser:


Step 1: Thoroughly wet the ends of the wrist straps where the small sleeves are situated. Use normal tap water. If the electric pulsing by the electrodes during use is too weak and not detectable, add small amounts of natural pure sea salt to the water or dip into colloidal silver (if you are user of colloidal silver) to increase connectivity with the skin.

Step 2: Find your pulse (main arteries on wrist, not the blue visible veins), position the straps with the sleeves between the straps and the skin. It is more comfortable to move around and perform normal tasks with the open end of the sleeve pointing in the direction of the elbow as can be seen in Fig. 1 above.


Position the sleeves around the wrists so that when the electrodes are pushed into the sleeves, the electrodes will be parallel and on top of your main arteries.

Push the electrodes into the sleeves as deep as it will go.

First method – position of electrodes on one arm:

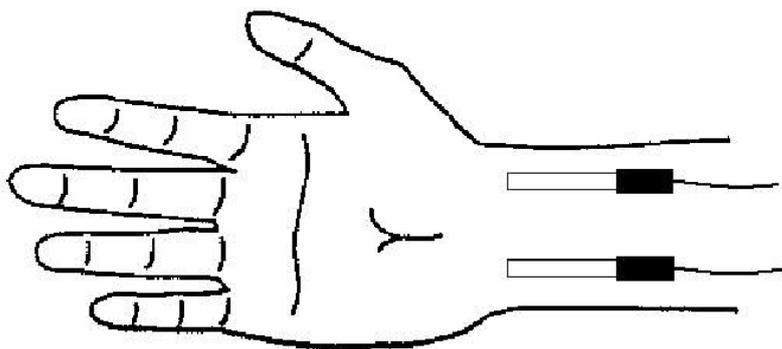

Second method – position of electrodes on both arms (main arteries and pulse easier to detect):

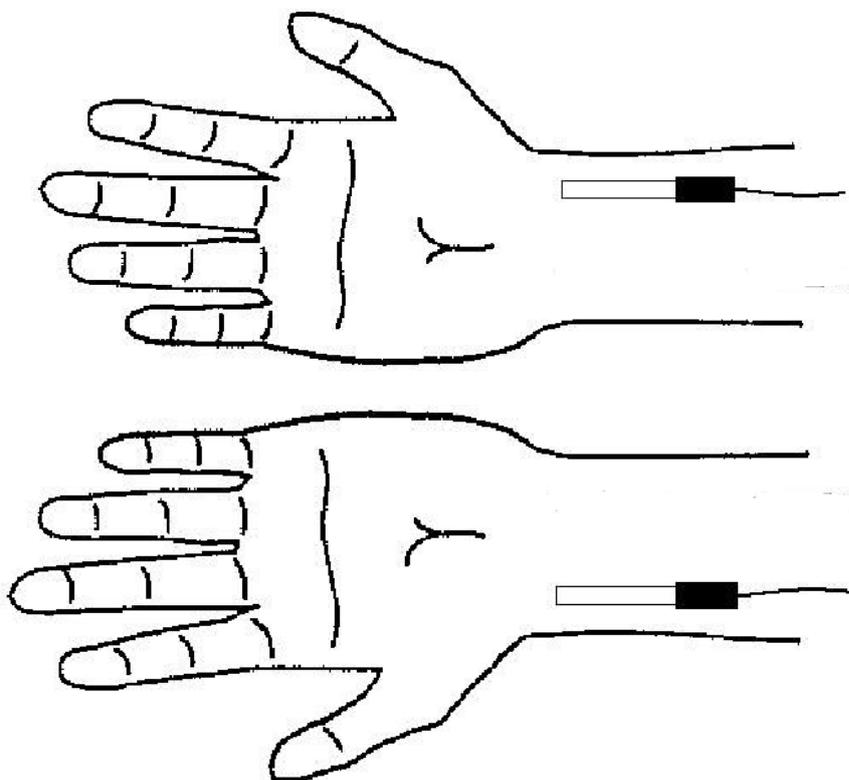

Step 3: Turn the amplitude (current intensity) potentiometer fully anti-clockwise to the minimum position

Step 4: Switch on the unit by pushing the small toggle switch on the front panel towards the potentiometer. The LED inside the transparent plastic box will light up continuously.

Step 5: Turn the potentiometer clockwise to approximately $\frac{3}{4}$ of rotation range until the LED begins to flicker. The LED flickers when the 3 - 7ma of switching current is flowing through your main arteries, now slowly turn the potentiometer further clockwise until you can feel the electric pulsing. Adjust to a comfortable level, not too high for comfort. It is not required that the pulse amplitude should cause a hurtful “stinging” sensation or twitching of the fingers.

Note:

- It is possible that some people may not feel, or barely feel the electric pulsing at full amplitude in the beginning of a session, your skin will become more conductive after a few minutes and you may need to adjust the intensity to a lower setting several times during a session.
- Whenever no current is flowing and the electrodes have lost contact, the LED will lit continuously even at maximum amplitude, check the electrode positioning.

Using the DC adaptor

The DC input is in the middle of the front panel. Any 3V to 24V DC (500mA min current) source with 2.5mm DC jack may be used. The supplied DC adaptor is normally 5V or 12V.

Plugging in the DC adaptor will light the LED regardless of the switch position. Move the toggle switch towards the “on” position (Towards the DC adaptor input to start the pulsing) for maximum current flow capability.

Batteries (non-chargeable, not supplied)

The unit has a current drain of below 50mA and can work for many hours using two 1.5v AA size batteries, Duracell Alkaline will offer the best performance. Alternatively two 1.5V rechargeable NiMh cells can be used. A 9V battery can be used (will not last as long as the 1.5V AA size batteries). Please note the polarity indication inside the battery compartment, **DO NOT INSERT BATTERIES WITH THE WRONG POLARITY** as it may damage the unit.

Rechargeable batteries (not supplied)

For recharging a 9V rechargeable battery (can be installed in battery compartment), plug in the DC adaptor, switch the toggle switch away from the DC input to the outside edge of the box.

A rechargeable NiMh 9V battery should normally be able to last for 1- 2 (1-hour) sessions.

Do not charge the battery if it has not fully discharged. The battery has run flat when no pulsing sensation can be detected anymore and/or the LED flicker rate becomes erratic.

Connect the DC adaptor to the unit. The LED will be on.

Switch the small toggle switch **away** from the DC adaptor input towards the edge of the front panel (indicated as “nu” in picture above).

Turn the amplitude potentiometer to the maximum clockwise position. The LED will flicker.

Short the electrodes (leaving them to touch each other). Turn the amplitude knob clockwise to maximum position.

Leave for a minimum of 12-15 hours to charge, do not leave on charge for more than a day.

Please note that rechargeable 9V batteries are not very powerful and may only provide enough power for one or two (1-hour) pulsing sessions. Rechargeable NiMh 1.5V (size AA) rechargeable batteries will offer much better performance, but should rather be recharged by taking the AA batteries out and charge them with a separate battery NiMH battery charger.

LED indication

Some users might be concerned that the LED is using too much power and that may drain the battery unnecessarily. This is not the case. The LED used is a selected Super-Bright LED that is utilised well below normal operating voltage (24V through a 82Kohm limiting resistor for 50% of the 4Hz pulse cycle) thus drawing on average only a negligible ~15 microamperes.

The Legal Stuff

The content provided in this publication is for informational purposes only and is not intended to diagnose or suggest a definite cure for any medical condition. Please consult your health practitioner if in doubt. Government regulators may say that testimonials are misleading and deceptive and not typical. If you believe the online testimonies to be honest, real life experiences, and would like to test this technology yourself, you may do so knowing that nothing is implied, promised or guaranteed.

Use this technology at own risk as a self-help natural healing method or as an experimental device for further research.

The manufacturer of the Eliminator Blood Pulsar and Colloidal Silver Generator is not a medical professional and offer these units solely as a replication of the Beck Protocol and all references to medical conditions and results are based on the research by self-help enthusiasts, self-application results and other publications by individual medical researchers.

For more information on the Beck Protocol and a database of testimonials regarding this technology: www.bobbeck.com

A simple explanation of how the protocol may work after studying many research materials:

The sub-atomic size electrons are like terminators on a mission, they are very focussed, they don't waste time and they move at light-speed from one end to the other following the shortest conductive route, they will not take a "detour" through the heart. Electrons don't get tired, trillions of them will pass through the blood each second, zapping viruses and bacteria along the way but not your body cells, which are electrically negatively charged and therefore electrons flow around them and not through them, so they are good friends with electrons at this limited current level and a low pulsing frequency of ~4Hz.

Viruses are criminals, invaders with a protein-like membrane or outer layer (fatty layer) and a spike receptor that is positively charged to attach to negatively charged body cells and immune cells, they inject their toxins into your body cells and will invade the body cell so they can multiply and feed on the nutrients meant to keep *you* healthy. The electrical field of the electron flow will change or destroy the protein layer of the virus membrane) and disables the bugger's ability to multiply or attach to body cells, hooray! The body's own

immune cells will get very excited and motivated by you helping them to disarm the enemy and will go into an intensified war effort on these de-weaponised criminals, they will kick them useless parasites out the door and you may experience some detoxing symptoms during the first few days or a week (slightly running nose, more visits to the loo to release the fluids, maybe even a slight headache for a day or two). You must drink lots of water, ozonised water is even better to help the liver and kidneys quickly expel the toxins from your body during this time. Start with a pulsing time of 10 minutes a day, listen to your body, gradually increasing the time over the span of two weeks until all detox symptoms and discomfort has left and you can now pulse for 1-2 hours.

What about bacteria? Well they are the easy prey and get barbecued on the spot! ☺.

Are you drinking Colloidal Silver? Now these heroes are negatively charged nano-sized particles that will attract the positively charged viruses to attach to them instead of your body cells, further contributing to the fight, a natural antibiotic, constantly increasing your immunity!

I personally suspect the silver colloids in the blood get further charged and excited by the electrical pulsing of the blood and becomes even more effective, it's really game over for the parasites! Keep silver colloid under the tongue for a minute before drinking a shot glass (25ml) once or twice a day (not the glass ☺) for faster absorption into the blood. Silver Colloid is reportedly absorbed into the blood before reaching the good microbes in the lower intestines, but if in doubt, take pro-biotic replenishments once in a while, no harm done.

While electro-pulsing at ~4Hz, the body cell membranes becomes permeable for a short while for nutrients to be absorbed and to expell toxins. During this period the body cells also invite the silver colloid comrades inside the cells for further house-cleaning and to expel toxins and parasites from your body cells. So, be careful if you are a smoker or consume those intoxicating alcoholic liquor during this time, it also gets amplified with increased absorption vcapability, why take toxins with the good stuff? Don't do it!

It will be better to not smoke, drink or eat shortly before or during pulsing (water is always welcome). Avoid cafeine intake shortly before or during pulsing as cell abosbtion will be amplified. 10 Minutes after pulsing you are of course welcome to eat and drink (preferably the healthy stuff) as normal.

STAY AWAY FROM GARLIC! Dr Beck reported that garlic might be poison to viruses and bacteria in the blood, but if absorbed into the body cells can be toxic to your cells as well.

A few personal experiences:

Friends have reported on the improvement of athritis, shoulder pains, hip pains, stomach pains, better clarity of mind and muscle pains relieve etc.

I have personally seen stomach pain or cramping disappear within 5 minutes of electro-pulsing of the blood on the wrists. That pestilence virus called Herpes causing sores/blisters appearing on the lips (or elsewhere): when you feel the tingling of that starting, zap that virus with one electrode (in a wet sleeve or wrapped in a wetted small piece of material) on the lip and on the spot where you feel the tingling, the other electrode (in a wet sleeve) below the lip for 10 minutes, three times a day. That blister will normally not develop further and his friend "The Cold" that normally takes advantage of your compromised immune system a week later, will get so scared when he see the on-spot execution of Herpes that he just don't stay around any longer, or maybe he got zapped because Herpes warned you in time to implement the blood pulsing. I also saw red swollen spots on the skin or face due to a

bacterial infection disappear (the swelling) within three minutes of localized zapping between the two electrodes. I experienced a discomforting pain below the stomach for three days while developing the unit, on day three I just got the circuit working on a test board, so I decide to put it to the test holding one electrode (in wet sleeve) on the spot at the front with the fingers and the other electrode positioned at the back against a cushion to keep it in place. I zapped for 45 minutes through the lower part of the body and the pain left within 45 minutes, it never returned to this day. I don't care what it was, it's gone!. I can't explain the speed of these healings in any biological way, but I thank Dr Beck for sharing his discovery!

Note:

You can make your own small sleeves or just wrap a soft strip of material around the electrodes in two to three layers and wet thoroughly for localised spot treatment.

When you experience healing, please be so kind to email your testimonial to:

pps69609@gmail.com

To your health and God bless!

Philip Smit

Indemnity Declaration

I, _____

ID No: _____

Hereby declare as follows:

1. I am aware that the manufacturer or supplier of the Blood Pulsing Device (also known as the Beck Protocol Blood Pulsar or The Eliminator Electro-pulsing Blood Cleanser) will accept no liability for any injury, adverse health related issues or death of the user directly or indirectly from the use of this device.
2. I understand the supplier of The Eliminator Electro-pulsing Blood Cleanser unit is not a medical professional and can as such not offer professional medical advice or any specific treatment recommendations.
3. I am aware that I should obtain medical advice from a medical professional regarding the use of an elctro-blood pulsing unit and the influence of such self-treatment on my current state of health.
4. I have done my own research regarding the Beck Protocol methods and devices online and have requested Mr Philip Smit, an electronic technician to construct the Blood Pulsing unit based on the original circuit published as freeware information by Dr Robert Beck, published on the internet with a few modifications to improve the user interface and to be able to also use the device connected to a DC adaptor connected to a 250V AC wall socket.
5. I am aware of the precautions regarding the use of the unit as published in the document that was downloaded here and specifically also referring to pages 10 and 11:
<https://www.bobbeck.com/pdfs/guide-to-using-the-beck-protocol.pdf>
6. I am using this unit as an experimental alternative or in addition to conventional medicine or treatment and understand that it is an experimental self-help device.
7. I am aware that there might not exist any medical clinical tests for this treatment.
8. I understand that the supplier or manufacturer of the device does not imply any guaranteed cure for any kind of illness whatsoever and I am using it purely as part of my own research and experimentation.

Signed at: _____

Signature : _____

Date Signed: _____

Please email signed form to pps69609@gmail.com